

8. Alarm Tone(s)



The Med-E-Lert Automatic Pill Dispenser features 3 different alarm tones and one silent setting to accommodate a number of different hearing levels. The display will show "A", "B", or "AB" depending on which tone is selected, and will show neither A nor B when silent mode is activated. To select an alarm tone, press and hold button "3" and simultaneously press button "1" to select a new tone. Continue pressing button "1" while holding button "3" to scroll through and listen to the different alarm tones.

▲ Press button 3 first, then press button 1 to choose the next alarm tone.

9. Testing the Med-E-Lert Functions ("CH" Mode)




The Med-E-Lert features a built-in testing function to ensure the device is working correctly. To begin test mode, press and hold button "1" for approximately 12 seconds until "CH" appears on the display. When the "CH" appears, the Med-E-Lert is ready to be tested. Press button "2" to perform a system function test. Once button "2" has been pressed, the Med-E-Lert will perform the following system tests:

▲ press button 1 for about 12 seconds, it display CH


■ CH 1 Display Test

The display will show "CH 1". All segments on the Med-E-Lert display will light up.

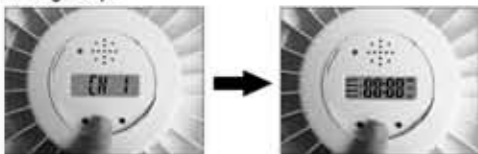


■ CH 2 Display Test

The display will show "CH 2". The motor will rotate the tray one compartment. If the tray does not rotate properly the motor will stop after a few seconds and the display will show the "CH 6" error message.

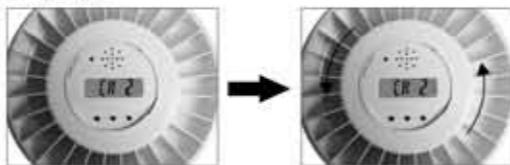

■ CH 3 Alarm Signal Test


The display will show "CH 3". The alarm will sound and the LED light will be lit. Turn the Med-E-Lert over to simulate dispensing pills and test the RESET function. Once the dispenser is turned over, the alarm and LED light should both be stopped. If the Med-E-Lert does not detect itself being turned over within ten seconds of starting the test, the display will show the "CH 4" error message.

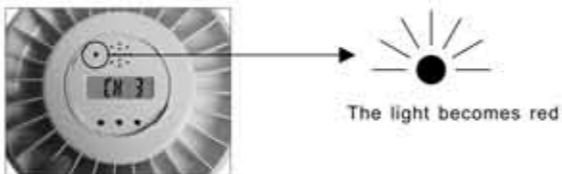

● To exit testing mode, press button "3". If there is a malfunction, remove and re-insert the batteries to reset the device.

10. Locking and Unlocking the Med-E-Lert

Insert the key into the device and turn it ¼ turn to the left to lock or ¼ to the right to unlock.

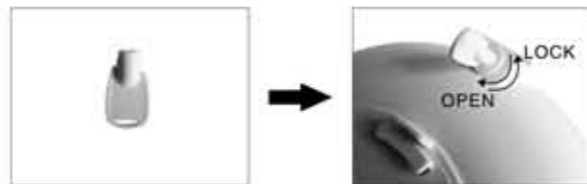

11. Stopping the Alarm



The alarm will stop once the Med-E-Lert senses that pills have been dispensed. Simply pick up the device and turn it over, carefully dumping the pills into your hand or a container.

12. Low Battery Alert



When the Med-E-Lert senses that batteries are low on power, an indicator will appear on the display (as shown). Replace the batteries immediately to ensure that the device continues to function properly.

Note: When the batteries are removed, the Med-E-Lert will save your alarm settings until new information is programmed. You do not need to reset your alarm(s) when batteries are replaced.



ActiveForever.com

Discount Therapy, Fitness, Pain Relief & Medical Equipment

1-800-377-8033

10799N. 90th Street
Scottsdale, AZ 85260

Warning: Do not place the MED-E-LERT in direct sunlight or strong ultraviolet light as it may cause malfunction in operation.

EASY SET-UP

MANUAL

1. Insert the Batteries



Insert the 4 "AA" Batteries into the compartment according to the guide imprinted in the bottom of the battery compartment and replace the cover, ensuring that it snaps into place.

2. How to Open the Medication Compartment



Press gently with your thumb on the front tab of Med-E-Lert as shown, then lift the cover.

3. Cleaning the Tray




Clean the medication tray using only a dry or damp cloth, cotton swab or sponge.


Note: Please make sure your hands are always as clean and dry as possible while using the Med-E-Lert in order to avoid the buildup of crumbs in the compartments or electronic components. Also avoid subjecting the Med-E-Lert to liquids or extreme temperatures.

4. Loading the Med-E-Lert Dispenser

Choose most appropriate disc, secure the disc in place, and fill the designated compartments.

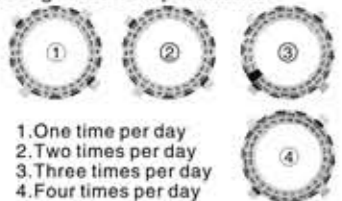

1. One time per day
2. Two times per day
3. Three times per day
4. Four times per day

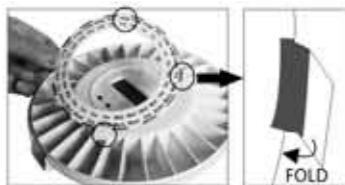

FOLD



The Med-E-Lert Automatic Pill Dispenser is supplied with 4 removable scheduling discs. The discs display the day of the week and the number of medication doses for each day (1, 2, 3, or 4 doses). Place the appropriate disc in the tray according to your desired number of doses per day, and fold the tabs down to secure the disc in the pill compartment. Use the day of the week and time of day marked on the disc as a guide for loading the dispenser.

Note: Assigning dosage times, organizing the medication in the compartments, and determining which medications go into each pill slot is the responsibility of you, your healthcare provider, or your caregiver.

The frequency of which the Med-E-Lert dispenser needs to be refilled depends on the number of medications taken per day. For example:

- One dose per day refill every 28 days
- Two doses per day refill every 14 days
- Three doses per day refill every 9 days
- Four doses per day refill every 7 days

5. 12/24-Hour Time Formats



▲ Press button 2 and button 3 simultaneously

The Med-E-Lert supports both 12 and 24 hour time formats, and automatically defaults to the 24 hour time format. To select the 12 hour time format, press the "2" and "3" buttons simultaneously and hold. The AM or PM symbol will appear when the clock is in the 12 hour format. Simply release the two buttons once you have selected the format you prefer.

6. Setting the Med-E-Lert Clock

(button 1- left; button 2- middle; button 3- right)



▲ Press button 1

1. Press button "1" and hold until the ":" symbol on the clock blinks, then release the button

BLINK



▲ Press button 2

2. Using button "2" for hours and "3" for minutes, set the correct time

BLINK
Set hour



▲ Press button 3

3. Press and release button "1" again to store the time

BLINK
Set minute

7. Setting and Checking the Alarm Times



1. To set the first alarm, press and release button "1". When the "Alarm 1" symbol will begin to blink, use buttons "2" and "3" to set the hours and minutes as you did setting the clock time.



2. Press and release button "1" again once you have selected the correct time to program that time into the Med-E-Lert. The "Alarm 2" symbol will begin to blink.



3. Repeat the above steps to program the remaining alarms. Once you choose the time for "Alarm 4", press button "1" again and the alarm symbols will stop blinking.

4. For less than 4 alarm times: if you wish to set less than 4 alarms, press and release button "1", and set the alarm time. After finished the last alarm you wish to set, then leave the display blinking. After one minute the Med-E-Lert will revert back to normal function



5. If you wish to cancel the alarm, press and release button "1", the alarm "1" symbol will begin to blink, and press button 2 until "--:" is displayed, which occurs after 23:00 or 11:00 p.m. then press and release button "1" to store the setting. To cancel alarm "2", alarm "3", alarm "4", using the same step.