




AsthmaMD™
Peak Flow Meter

Empower your asthma patients.

Teleflex is improving respiratory care with every breath.

Inspired by physicians and clinicians who help develop products, we are proud to offer the portable, pocket-sized AsthmaMD™ Peak Flow Meter to help patients measure their lung performance and anticipate attacks.

- Measures peak expiratory flow
- Lightweight, portable, and accurate
- Can be easily cleaned without disassembly
- Compatible with a convenient physician-developed mobile app

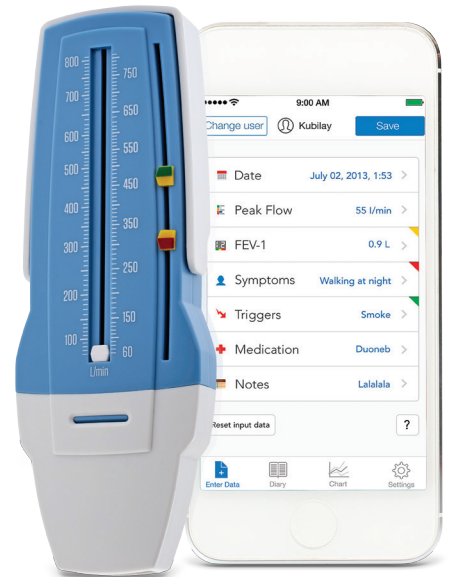

Benefits




The Clinician
Accurate, precise,
and designed for
children and adults¹



Your Institution
Meets all ATS standards
for portable Peak
Flow Meters²



The Patient
Durable, lightweight,
and easy to use

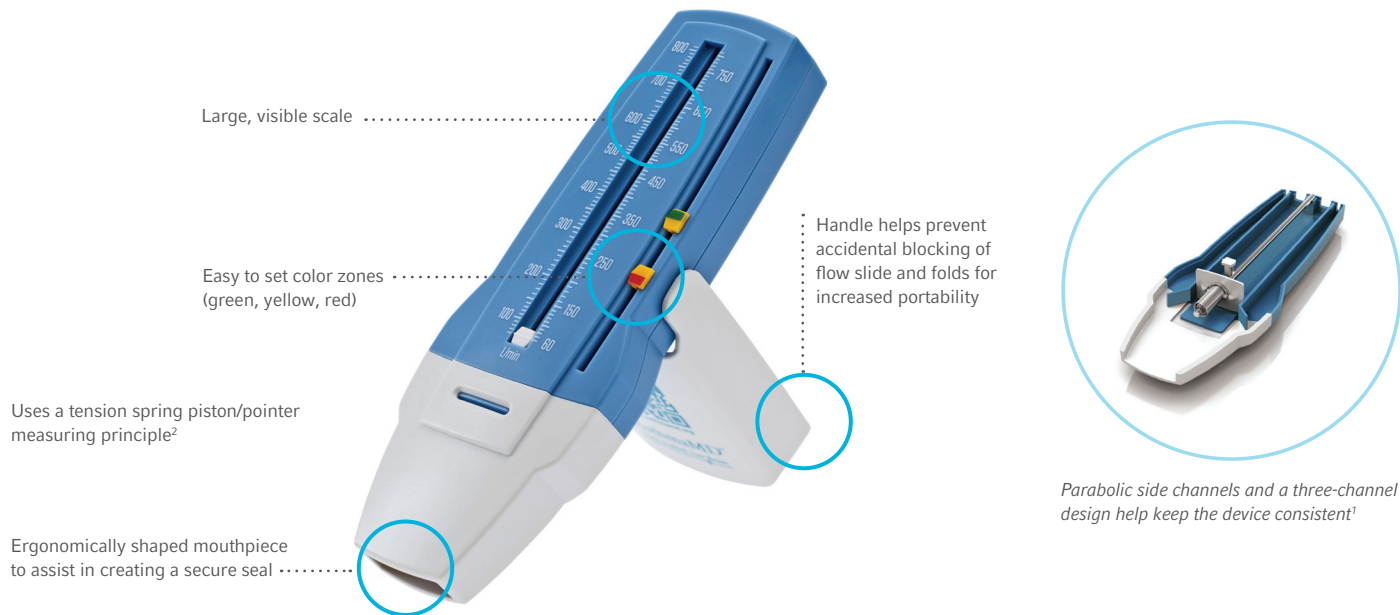


Parabolic side channels and a three-channel design help keep the device consistent¹

AsthmaMD™ Peak Flow Meter Fast Facts

- Compatible with a convenient app for tracking peak expiratory flow*
- Universal 60–810 LPM range
- Colorful, integrated sliding zone indicators
- Does not contain natural rubber latex or DEHP¹
- Meets all ATS standards for portable Peak Flow Meters²

AsthmaMD™ Peak Flow Meter Mobile App*

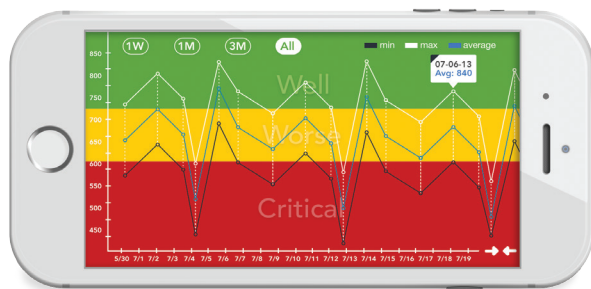

The AsthmaMD™ Peak Flow Meter companion app is an easy-to-use app that logs peak flow meter readings, tracks triggers, symptoms, medications, and notes. This allows the patient to replace paper journals and create a digital overview of their asthma activity.*



*DISCLAIMER: The AsthmaMD™ mobile app is a third-party app offered by AsthmaMD, Inc. Teleflex has no control over such third-party app, nor does Teleflex review, approve, monitor, endorse, or make any representations or warranties with respect to such third-party app. Phone also not included.

AsthmaMD™ Peak Flow Meter

ITEM CODE	DESCRIPTION	CASE QUANTITY
1802	AsthmaMD™ Peak Flow Meter	10

Specifications

Range	60–810 l/m (ISO23747, ATS)
Scale Increment	10 l/m
Accuracy	<+/- 10% or 10 l/m
Design	Single-Patient Use
Standard Compliance	ISO 23747, ATS 26 Flow-Time Waveform Testing, FDA 820
Cleaning Instructions	Clean at least once every two weeks. Hand wash: wash with warm water and mild liquid soap. Rinse thoroughly and air dry completely prior to use. Do not boil AsthmaMD™ Peak Flow meter. With proper care the AsthmaMD™ Peak Flow meter can last up to 2 years ¹

To place an order call 1.800-397-5899



[Click for price and product details](#)

References

1. Data on File, Genia Medical, Inc. 2. 510(k): K101380

AsthmaMD is a trademark of Genia Medical, Inc.

Teleflex, the Teleflex logo, and Hudson RCI are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this document may not be available in all countries. Please contact your local representative. All data current at time of printing (01/2017). Subject to technical changes without further notice.

© 2017 Teleflex Incorporated. All rights reserved. MC-002358